


SERVOPEN®

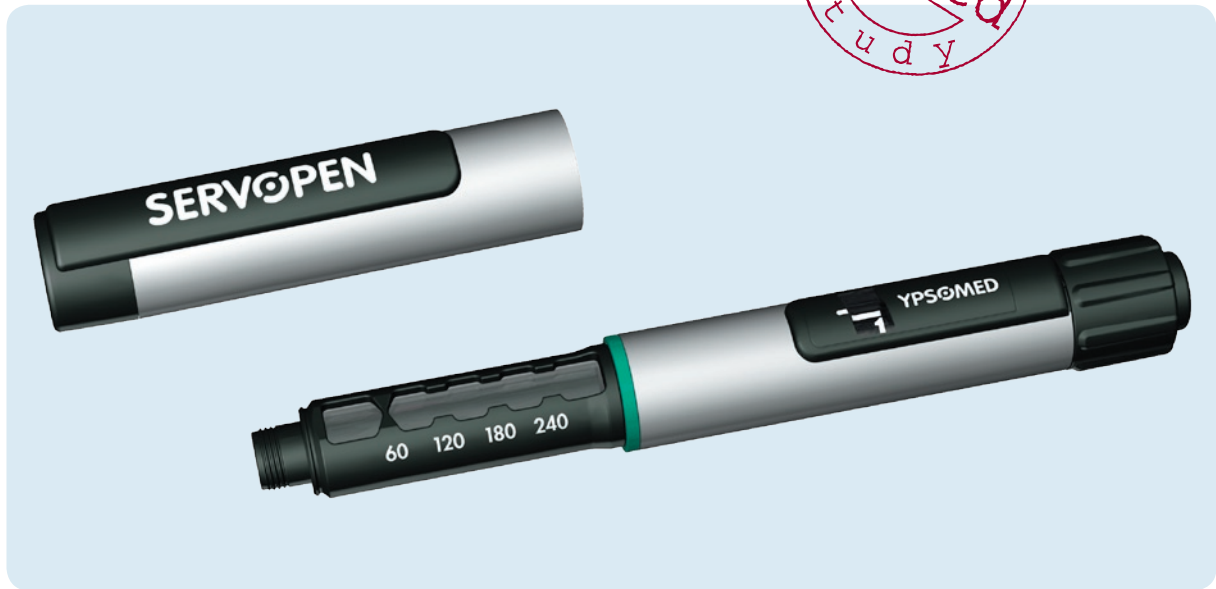

Automatic Variable Dose Pen for Insulin and Other Therapies


High quality reusable pen that sets new standards in terms of ease of use and patient friendliness. Through the spring-assisted delivery it offers the shortest available stroke distance on the market. Time required for cartridge exchange is minimised.

Multidose	Variable dose
Monodose	Fixed dose


Features and benefits

- Spring-assisted delivery: Low force required, shortest available stroke distance on the market
- Real-time display: Possibility to interrupt injection and to know how much still has to be delivered
- Very large numbers: Good readability even for visually impaired patients
- Simple cartridge exchange: Push-back piston rod and simple bayonet connection for cartridge holder
- Material: High quality look & feel including metal outer materials
- Dosing: Variable dosing in single steps from 1–60 IU, with dose-correction function
- Audible feedback during dose dialling and injection
- Stop-last-dose (stop-300): The dialled dose cannot exceed the volume remaining in the cartridge
- Exchangeable colour rings: To colour-code pens when using different insulins


Main handling steps

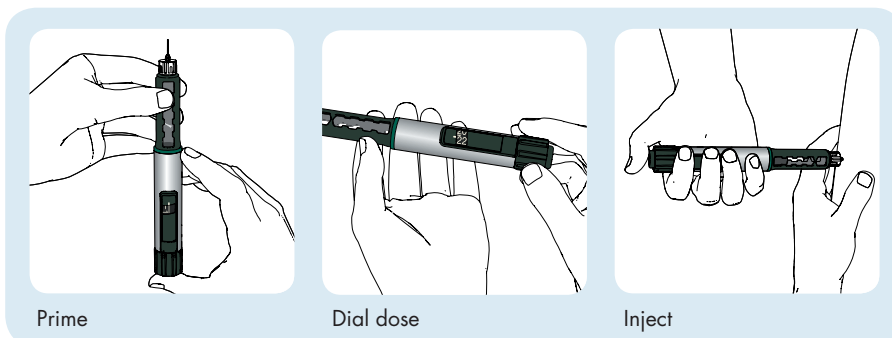

Pen dimensions

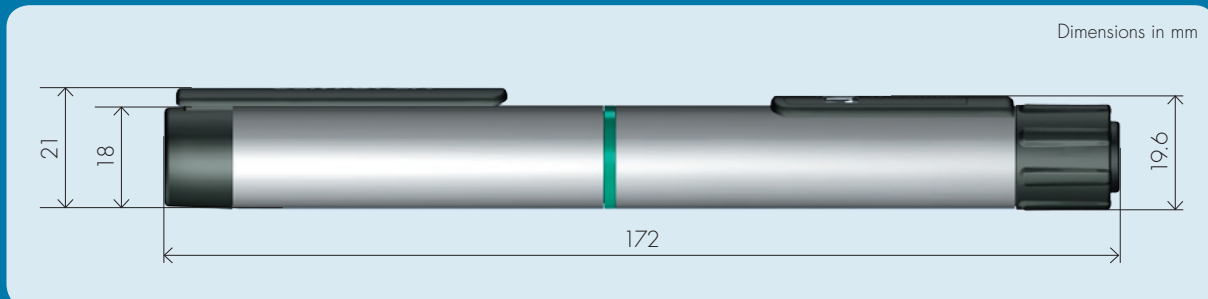

Cartridge type

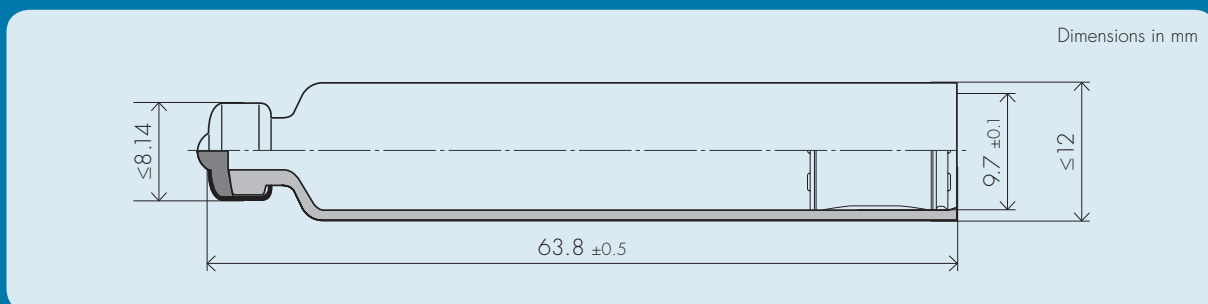

Standard 3ml cartridge

COMPATIBILITY BETWEEN EACH NEW CARTRIDGE AND THE SERVOPEN® REQUIRES VERIFICATION

Pen specifications

- Min. dose: 1 IU (10 μ l)
- Max. dose: 60 IU (600 μ l)
- Number of dose increments: 60
- Compatible with standard and safety pen needles – Ypsomed Penfine®/Clickfine® and Clickfine® AutoProtect™ pen needles with unique click-on function are recommended and can be prepacked with every pen.