

# UnoPen™ Fix



## Fixed Dose Disposable Pen


Delivery of multiple fixed doses made simple with the proven pull-push dose delivery concept. Clear audible and tactile click signals reassure patients and ensure correct handling. State-of-the art gearing enhances handling comfort during dose delivery.

Multidose	Variable dose
Monodose	Fixed dose

## Features and benefits

- Simple and easy handling with pull – push principle
- Clear visual and audible signals confirm completion of “dose setting” and “injection”
- Geared mechanics enhance handling comfort during injection
- Can be customized for a wide range of cartridge sizes and types

## Main handling steps

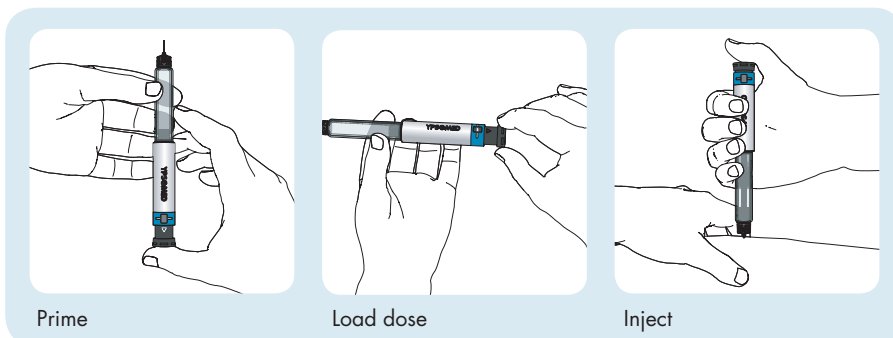


Prime


Load dose

Inject

## Pen dimensions

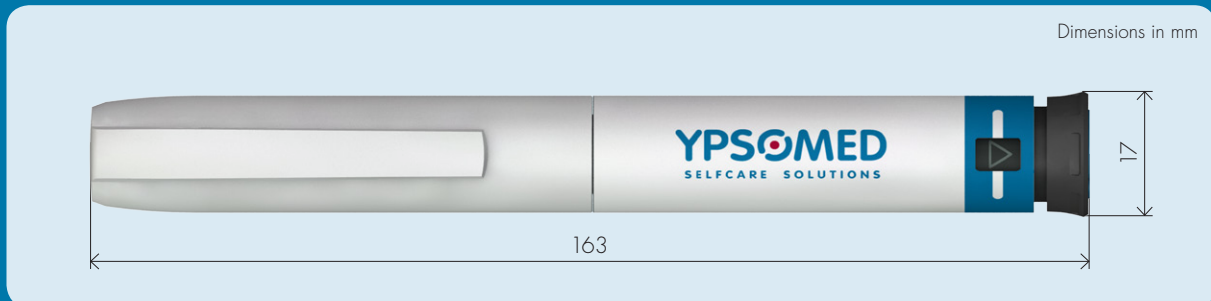

## Cartridge type

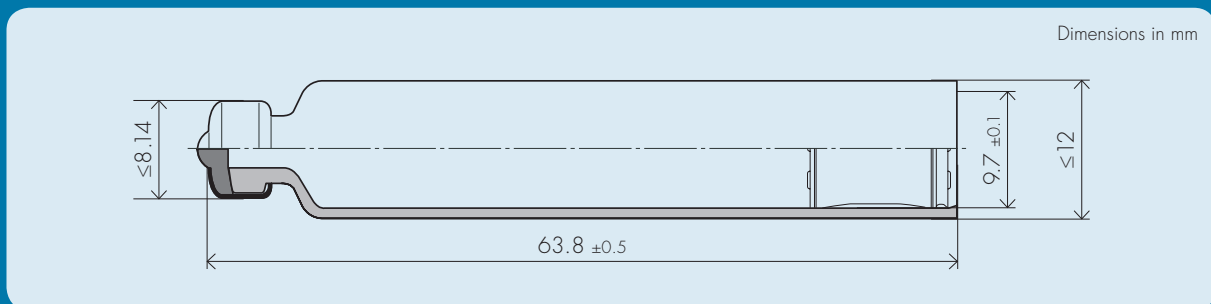

Standard 3ml cartridge

COMPATIBILITY BETWEEN EACH NEW CARTRIDGE AND THE UnoPen™ Fix REQUIRES VERIFICATION

## Pen specifications

Note: The maximum dose impacts the device length. The shown pen is designed for a dose of 0.1 ml

- Dosing range: 0.04–0.2 ml
- Compatible with standard and safety pen needles – Ypsomed Penfine®, Clickfine® and Clickfine® AutoProtect™ pen needles with unique click-on function are recommended and can be prepacked with every pen.