


# UnoPen™ Memo

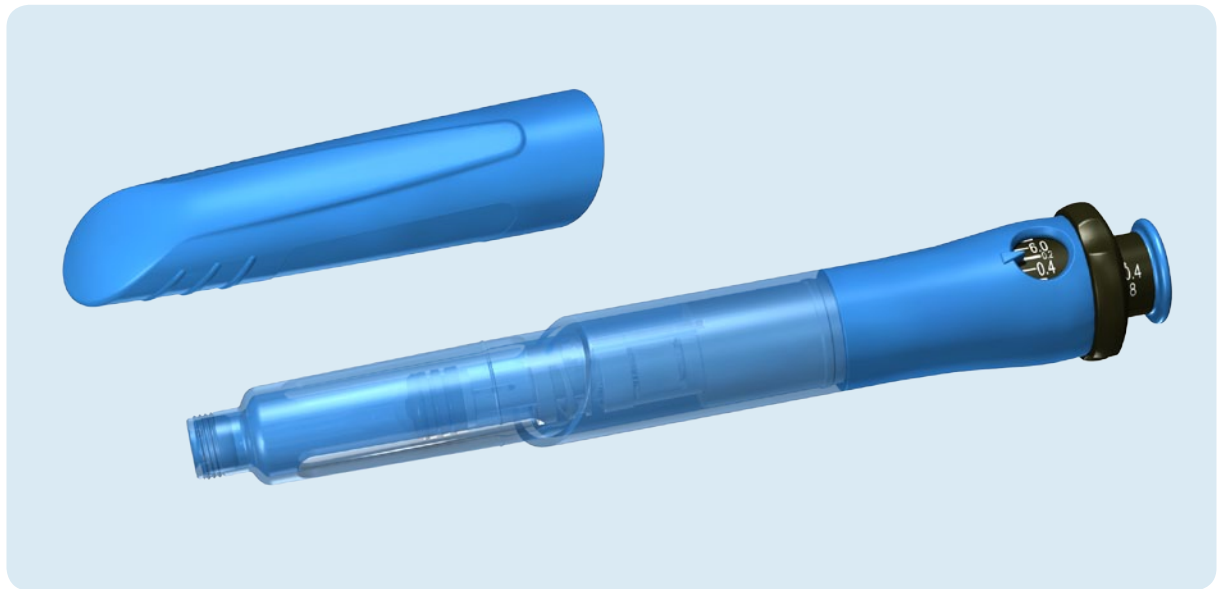

## Variable Dose Disposable Pen with Dose Memory


Variable multi-dose pen with unique dose memory feature. A benefit for patients who need to repeatedly inject the same dose over a prolonged time period. Patients only need to set the dose once for each new cartridge. This simplifies handling and reduces the risk of dosing errors. The pen gearing enhances handling comfort during loading the dose injecting.

Multidose	Variable dose
Monodose	Fixed dose


### Features and benefits

- Dose memory function for repeated injections of the same dose
- Easy and intuitive to operate by simply loading the dose and injecting
- Geared mechanics enhance handling comfort during injection
- Simple dose selection (up to a maximum of 30 doses possible)
- Can be customized for a wide range of cartridge sizes and types, including dual chamber cartridges


### Main handling steps

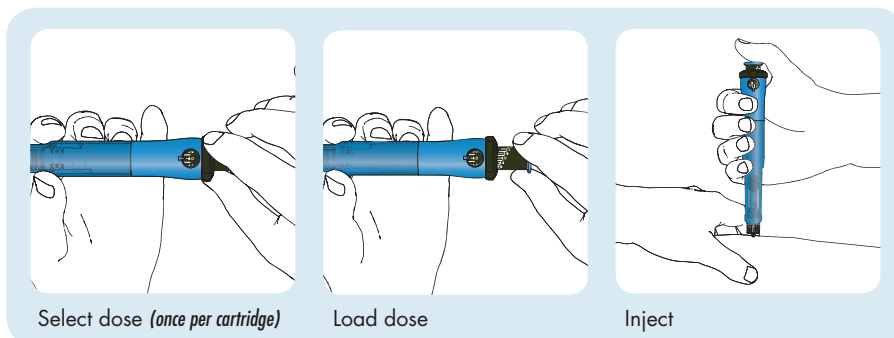

## Pen dimensions

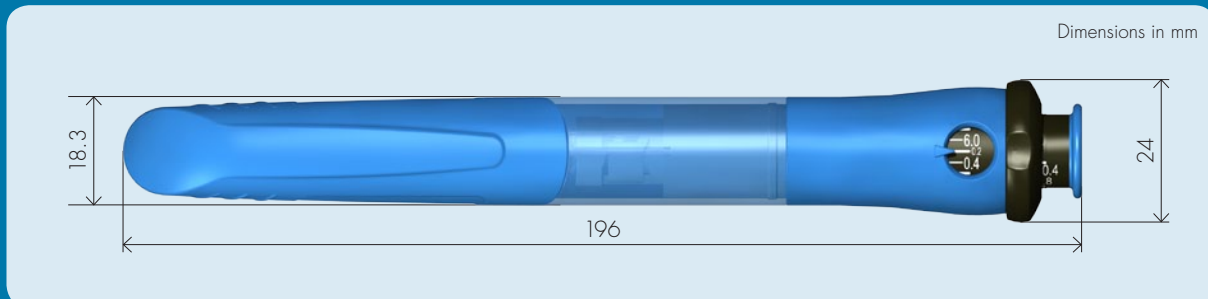

## Cartridge type

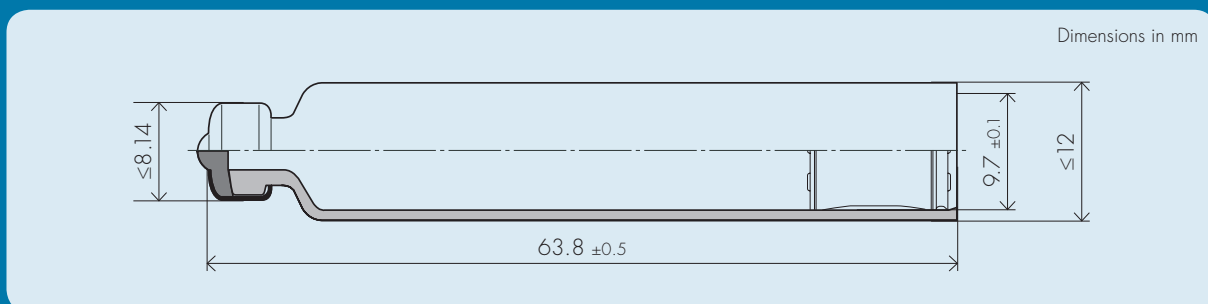

Standard 3ml cartridge

COMPATIBILITY BETWEEN EACH NEW CARTRIDGE AND THE UnoPen™ Memo REQUIRES VERIFICATION

## Pen specifications

Note: The maximum dose impacts the device length. The shown pen is designed for a dose of 0.3 ml

- Dosing range: 0.01–0.6 ml
- Max. number of dose steps: 30
- Compatible with standard and safety pen needles – Ypsomed Penfine®, Clickfine® and Clickfine® AutoProtect™ pen needles with unique click-on function are recommended and can be prepacked with every pen.